




AKJA 2000 DIVISIBLE

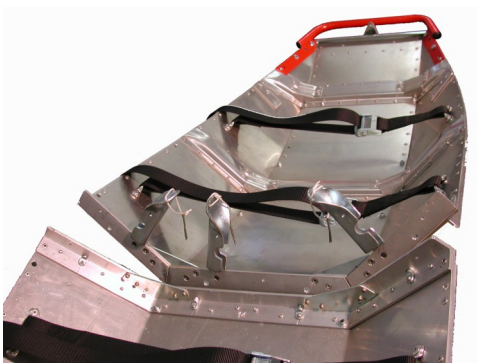

The classic Akja - Rescue Sled in **divisible version** in length of 200 cm.

Designed for **2-rescuer** operation in **all kind of snow** covered terrain, groomed or ungroomed.

Divisible in two halves for transport on the back in unexploited areas.

Long lasting durability: Manufactured completely from high-strength, corrosion-resistant aluminum alloy for many years of reliable, maintenance-free action.

Solid and **quick dividing - mechanism** provides a safe hold while riding and quick parting in the middle if necessary for transport.

Stable running characteristics: Reinforced steel runners with especially hardened fins provide stable running even in very icy conditions.

Two pairs of poles (for two rescuers), can be detached easily by bayonett-lockings. The grips of the handles are high enough for tall rescuers. With integrated stabilization profiles for increased bending strength.

Optionally, divisible handles are available.

Bayonett-lockings are wear-free, 100% icing-safe and easily to handle.

Carrying grips on both ends of the sledge for comfortable lifting and carrying of the sled.

Four fixation belts with aluminum buckles for safe fixation of the injured; resistant against breaking and icing, can be operated with skiing-gloves.

Art. No.	93032	
Dimensions	Length:	2000 mm
	Width:	560 mm
	Length divided:	1.000 mm
Weight	Total incl. poles:	27 kg
Comments	Divisible handles optionally available.	