




## AKJA 2000—SHORT VERSION

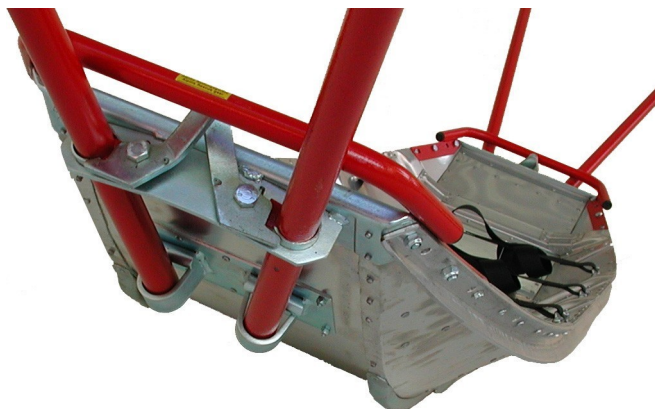

Designed for **2-rescuer** operation in **all kind of snow** - covered terrain, groomed or ungroomed.

The **classic Akja - Rescue Sled** in one-piece version. Time-proofed design and workmanship, thousands of them worldwide running.

Especially good in **steep and/or icy slopes** and moguls.

Due to the **shorter length** this model is very practicable for off-piste rescue operations.

Also sometimes for **transport with a gondola** this shorter version is needed.



**Long lasting durability:** Manufactured completely from high-strength, corrosion-resistant aluminum alloy for many years of reliable, maintenance-free action.

**Stable running characteristics:** Reinforced steel runners with especially hardened fins provide stable running even in very icy conditions.

**Two pairs of poles** (for two rescuers), can be detached easily by bayonett-lockings. The grips of the handles are high enough for tall rescuers. With integrated stabilization profiles for increased bending strength.

**Bayonett-lockings** are wear-free, 100% icing-safe and easily to handle.

**Carrying grips** on both ends of the sled for comfortable lifting and carrying of the sled.

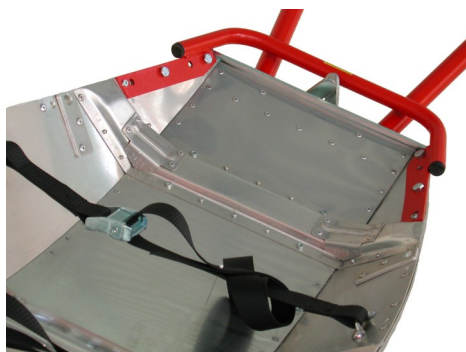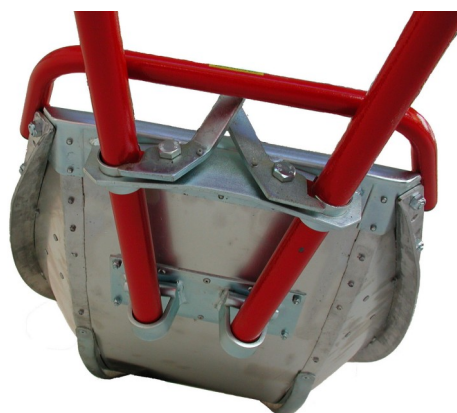

Art. No.	<b>93059</b>	
Dimensions	Length:	2.000 mm
	Width:	560 mm
	Length of poles:	1.650 mm
Weight	Total incl. poles:	26,0 kg
Comments	Short version of Art. 93031	

