



## ADJUSTABLE SUSPENSION TYRAL




Loops leg area adjustable from 30cm to 150 cm. Black 21mm  
Loop head area fixed length of 100 cm. Yellow 21mm

Attachment to the stretcher using Twistlock carabiners.  
The suspension on the carabiner is hung on the rope using a  
rigger plate.

The suspension can be converted from a horizontal to a  
vertical suspension without tools.  
Simply unthread the two yellow ribbon loops from the  
OpenRing and use them for vertical lowering.

All parts certificated: EN 354, EN 566 und EN 795-B.

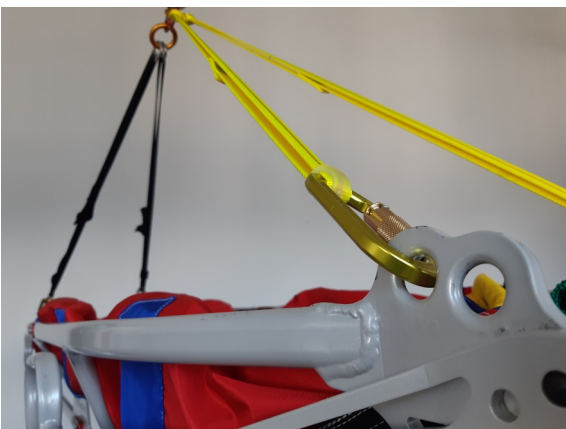

Art. No.	<b>93258</b>
Dimensions	Length Head: 100 mm Length Feetarea: 30-150 mm
Weight	Total: 1,4 kg
Comments	Delivered with TYRbag