



# GLIDING-SKI-MODUL TYRAL



The two Gliding-Ski-Modules are also **easy to apply** in soft snow and allow the TYRAL stretcher to remain easy to slide and pull even in muddy snow.

The identical and symmetrical shape avoids confusing the attachment.



The Gliding-Ski-Modules can be attached in **three different** positions.

- For flat terrain, apply offset to the right and left (head and foot section).
- Attach to both sides of the foot for vertical downwards.
- Attach to both sides of the head for vertical upwards.

The two sliding ski modules can be easily attached to the stretcher **without tools** using a pipe clip.

Another advantage is that the two Glide-Ski-Modules can remain on the stretcher even when the wheel is in use - no constant conversion necessary. The wheel support is only removed for crossing snow fields.

A plastic tub was deliberately avoided in order to be able to offer a **sustainable, long-lasting** product.

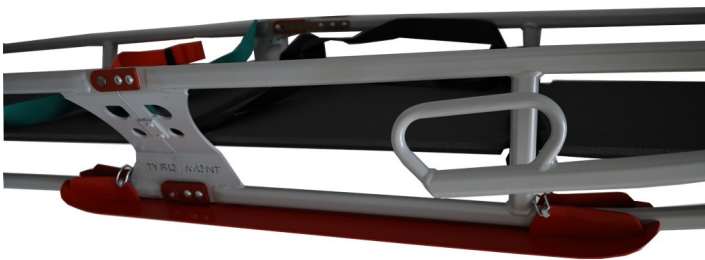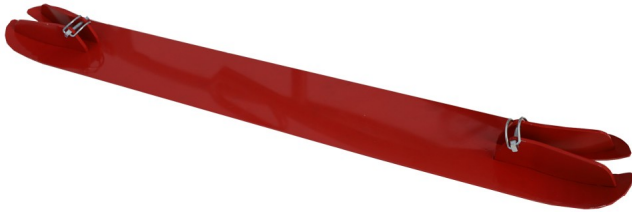

Art. No.	<b>93262</b>	
Dimensions	Length:	935 mm
	Width:	120 mm
Weight	Total per pair:	3,3 kg
Comments		

